



# VS 200

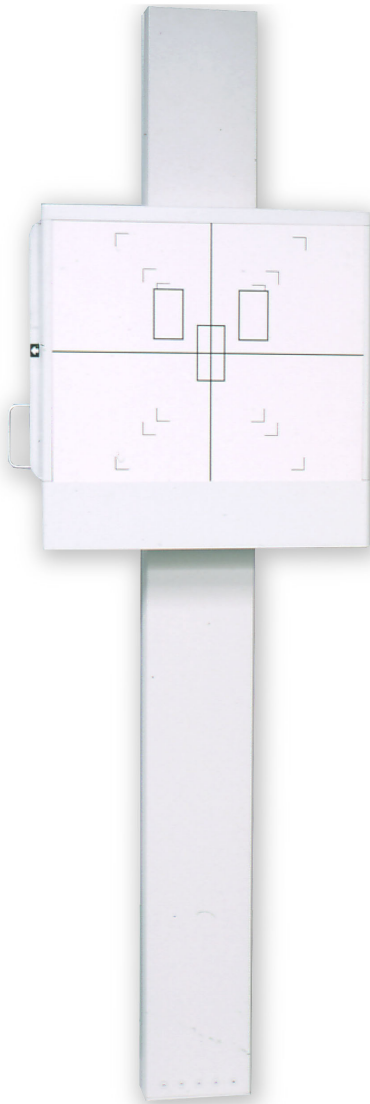

.....Vertical Wallstand.....

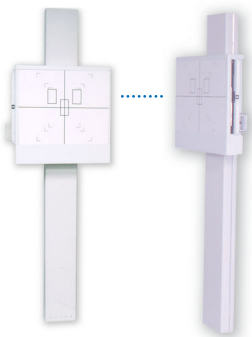

## VS 200 Vertical Wallstand

Designed specifically for clinical upright radiography, the VS 200 Vertical Wallstand is a stable film platform. The VS 200 meets all the requirements for upright radiography. The image receptor lowers to floor level to permit weight bearing studies.

The slender column design, structured architectural look, and the concealment of the cabling within flexible tubing all contribute to the VS 200's professional appearance. The VS 200 allows sturdy, smooth, vibration-free movement and provides long-lasting quality and durability.

Ideal for use in a hospital, clinic or private practice, this wallstand features a simple, spring-loaded mechanical lock for vertical travel. Minimal in-field installation required.

### System Options:

- ▲ True-speed buckys
- ▲ Trayless autoload buckys
- ▲ Del Medical 17" x 17" grid cabinets
- ▲ Automatic exposure control
- ▲ Positive beam limitation
- ▲ Patient hand grips (P/A)
- ▲ Patient hand grips (overhead lateral)
- ▲ Full line of high quality grids

### Technical Specifications

Height (including vertical travel):	84" (213cm)
Depth:	12 1/4" (31cm)
Width:	25" (64cm)
Weight (with image receptor):	200 lbs (91kg)
Front panel aluminum equivalency:	less than 0.7mm

### Comprehensive Equipment Warranties

Complies with all CDRH performance standards in effect at the date of manufacture.

Specifications subject to change without notice.

Made in U.S.A.

UL-Listed

